



SCINTILLA®  
CARE

# COMFORT CARE BEDS



# CARE BEDS

**SCINTILLA<sup>®</sup>**  
**CARE**

Are you looking for a care bed that offers both stability and a variety of functions? We have exactly what you need! Our care beds are the perfect solution for your needs.

With its robust construction and high-quality materials, our care beds offer exceptional stability and safety. In addition, our care beds are equipped with a variety of functions that make comfort and care easier.

Contact us today for more information!

[info@funke-medical.de](mailto:info@funke-medical.de)

High level of safety...  
...unlimited comfort



# Convincing facts and features

---





**Modern and extremely stable:**

Our care beds are available in two versions, the **Scintilla 1000** and the **Scintilla 2000**.

Both models are characterised by their enormous stability, which is guaranteed even in the highest position.


**Safety and protection:**

Both models can be fitted with full-length or separable side rails.

The full-length side rails can be raised to the highest position with little effort. The separated versions have an intermediate stop, which can be used as a stand-up aid or as a guide for the patient. The side rails are available in three versions:



Full-length side rails

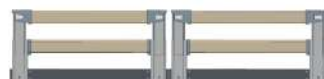

Separable side rails  $\frac{1}{2} + \frac{1}{2}$



Separable side rails  $\frac{2}{3} + \frac{1}{3}$





Incomparably solid...



# Scintilla 1000

## **Comfort:**

The height of the care bed can be adjusted thanks to the very stable X-scissor design. Even in the highest position the care bed ensures stability. The care bed can be easily moved to any position using the remote control. The **Scintilla 1000** provides support when getting up and helps with the transfer. To offer even more comfort, the grid frame is also available with an optional ergonomic panel system.

## **Work ergonomics:**

The **Scintilla 1000** care bed is a safe and stable bed that can be adjusted conveniently to the correct working height. The remote control can be locked with a key. The backrest has a virtual pivot point that prevents shearing forces on the coccyx. All additional options as shown in the catalogue can be attached to the care bed on request. All Funke Medical care beds are also double braked. This makes it easy to move a bed that is tilted in the room.

## **Economic efficiency:**

The **Scintilla 1000** is made from high-quality materials and guarantees a long service life. The adjustment options of the **Scintilla 1000** are designed to increase the patient's independence. The care bed can be customised to suit current needs. The **Scintilla 1000** fulfils the European legal regulations and is CE certified according to the MDR 2017/745. In addition, the care beds are TÜV tested. Our care beds do not only offer comfort and stability, but can be customised to suit any room and taste due to our various design options.





# Scintilla 1000 - technical specifications



Technical specifications	Scintilla 1000
External dimensions bed frame	90 x 200 cm
External dimensions	103 x 207 cm
Dimensions lying surface	90 x 200 cm
Bed extension	10 cm at head end
Optional bed extension	10 cm at foot end
Height adjustment	38 - 80 cm
Intermediate stop	Not available
Adjustment of back section	70°
Adjustment of thigh section	30°
Adjustment of lower leg section	22°
Trendelenburg	Not available
Anti-Trendelenburg	Not available
Brakes	Standard: Axle braking system (2 by 2)  Optional: Central braking systems (4 wheels)
Max. load	225 kg
Max. patient weight	175 kg
Net weight of the bed	115 kg
Optional full-length side rails	yes
Optional separable side rails	yes
Accessible to the lift	yes
Optional wheel covers	no
Possibility to fit additional options	yes
Electrical connection	100 - 230 V AC

# Function... ...meets design



Piemont Mokka

Country Oak

Beech

## Dimensions:

The **Scintilla 1000** care bed is available in the size 90/200 cm.

## Bed extension:

The **Scintilla 1000** is 2 meters long as standard and can be extended by 10 cm at both the head end and the foot end. The separated side rail can be extended with an intermediate piece.


The **Scintilla 1000** care bed is available in three wood colours.



Your care...  
...is our mission



# Scintilla 2000 - perfect symbiosis of function und design

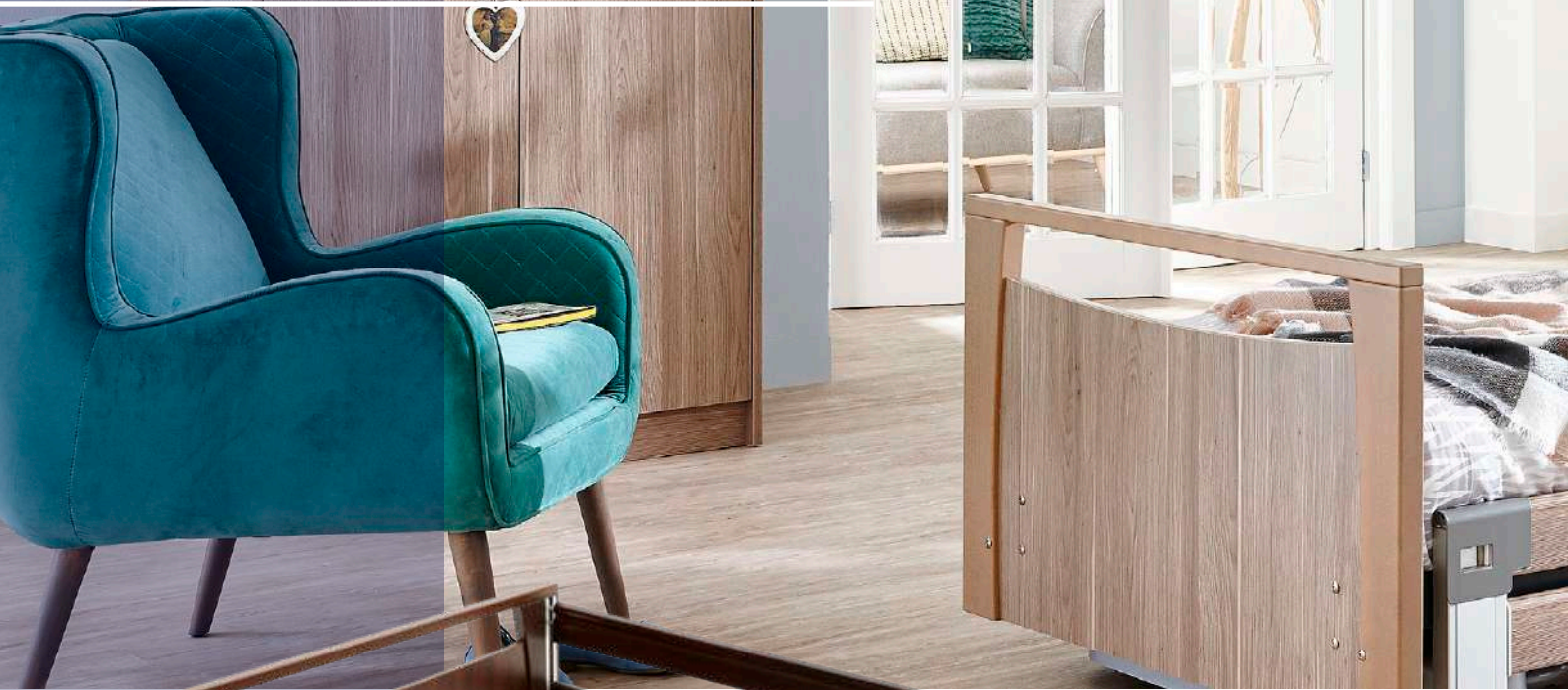



### Comfort:

The **Scintilla 2000** care bed can be adjusted using a very stable Y-scissor construction. Thanks to a technical refinement, the care bed can be transformed into a luxurious armchair. In this comfort setting patients can, for example, have breakfast or watch television in an upright position. The care bed also offers a low access height of only 25 cm and the luxurious cover on the wheels creates an elegant flair.

### Work ergonomics:

With its ease of use and many functional options, the **Scintilla 2000** makes care easier. Special attention has been paid to patient safety with this care bed. In addition, this care bed can be equipped with a 10-button remote control that allows the anti-Trendelenburg function and the Trendelenburg function to be set using a single button. Of course, additional options are available for this care bed.

### Economic efficiency:


The high-quality materials of the **Scintilla 2000** guarantee a long service life. Thanks to a variety of options the **Scintilla 2000** care bed makes it possible to take care of a larger group of patients. The **Scintilla 2000** fulfils the European legal regulations and is CE certified according to the MDR 2017/745. In addition, the care beds are TÜV tested. Besides safety and stability, the **Scintilla 2000** has a design that blends harmoniously into the room.

# Scintilla 2000 - technical specifications

Technical specifications	Scintilla 2000
External dimensions bed frame	90 x 200 cm XL: 110 x 200 cm
External dimensions	103 x 207
Dimensions lying surface	90 x 200 cm
Bed extension	10 cm at head end
Optional bed extension	10 cm at foot end
Height adjustment	25 - 80 cm
Intermediate stop	Standard
Adjustment of back section	70°
Adjustment of thigh section	30°
Adjustment of lower leg section	22°
Trendelenburg	14° locked
Anti-Trendelenburg	14°
Brakes	Standard: Axle braking system (2 by 2)  Optional: Central braking systems (4 wheels)
Max. load	225 kg (XL: 325 kg)
Max. patient weight	175 kg (XL: 290 kg)
Net weight of the bed	125 kg
Optional full-length side rails	yes
Optional separable side rails	yes
Accessible to the lift	yes
Optional wheel covers	yes
Possibility to fit additional options	yes
Electrical connection	100 - 230 V AC



# Unmistakably... ...versatile

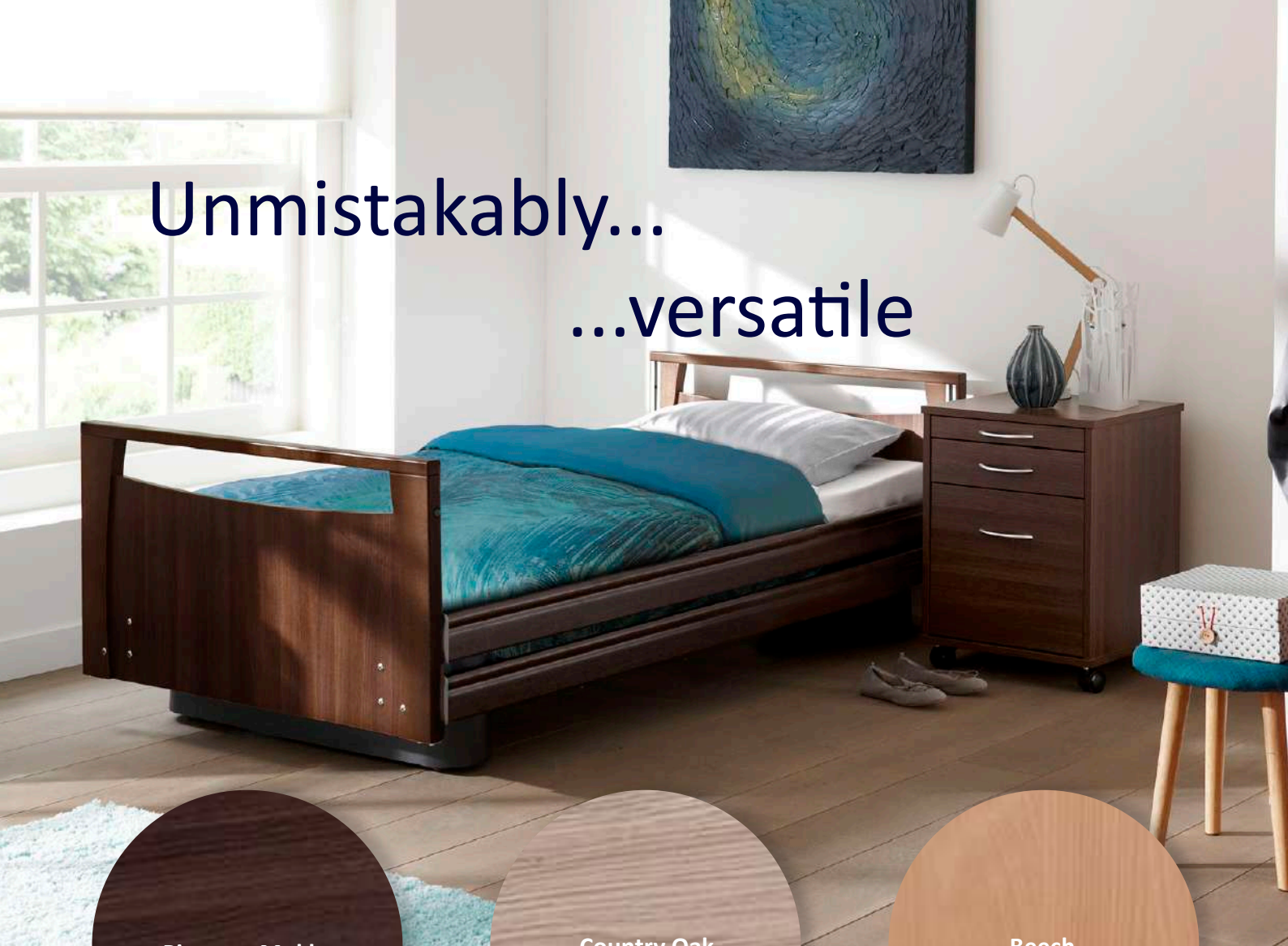

Piemont Mokka

Country Oak

Beech

## Dimensions:

The **Scintilla 2000** is available in the size 90/200 cm. Additional to the regular size, the care bed is also available in the oversized 100/200 cm version.

## Bed extension:

The **Scintilla 2000** is 2 meters long as a standard and can be extended by 10 cm at both the head end and the foot end. The separated side rail can be extended with an intermediate piece.

The **Scintilla 2000** care bed is available in three wood colours.

# Accessory

---



**Bedside lamp**



**Automatically/  
manually adjustable  
trapeze bar**



**Comfort lying surface  
with panel system**





**Stand-up bracket**




**Underbed lighting**

**Further options:**

- Battery
- Wall spacer
- Infusion stand



# Furniture

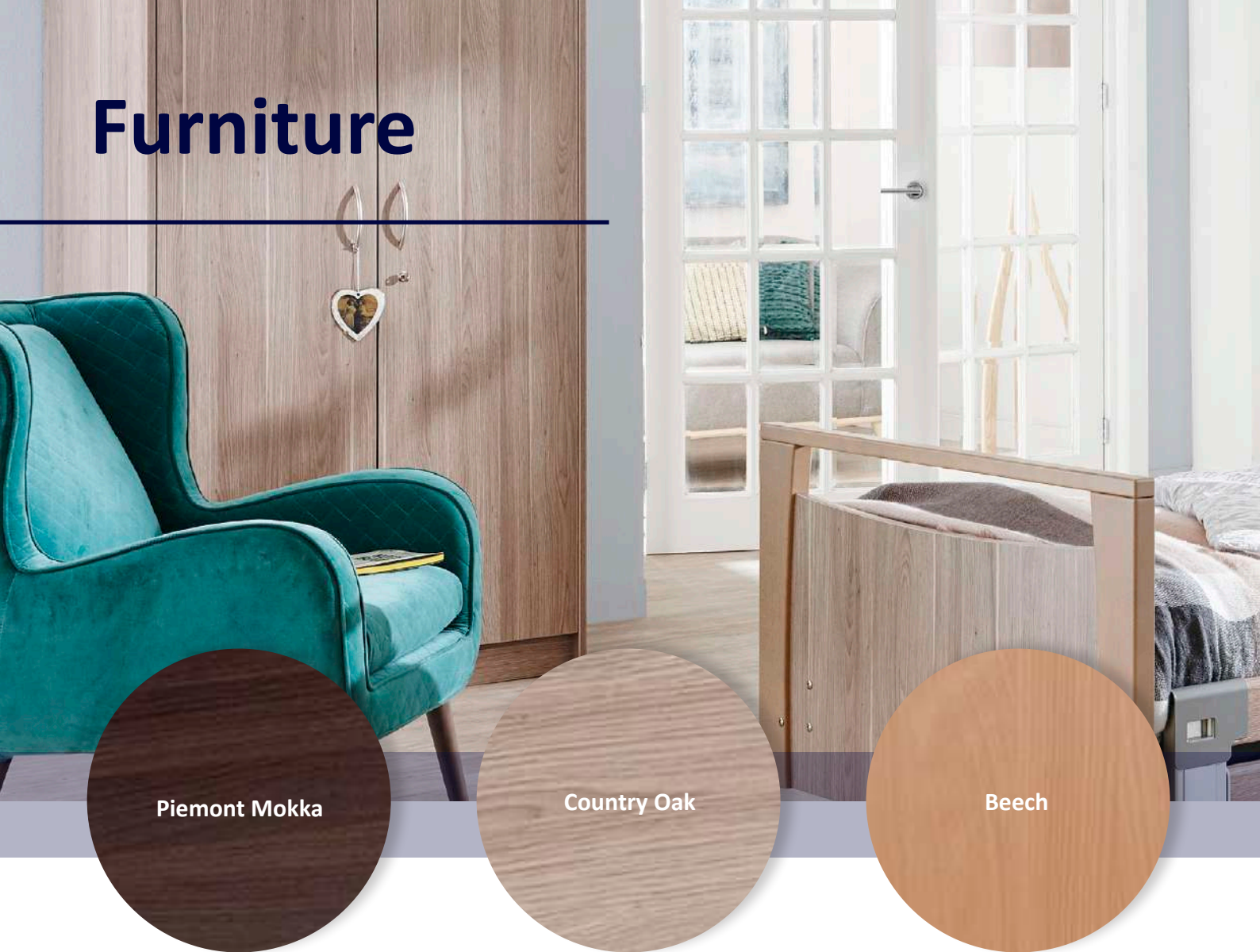

Piemont Mokka

Country Oak

Beech

To complete the interior, you can choose matching bedside tables and wardrobes from three different wood colours.

## Bedside table 1:

- Dimensions 48 x 48 x 74 cm (W x D x H)
- 4 wheels



## Bedside table 2:

- Dimensions 53 x 46 x 87 cm (W x D x H)
- Hydraulically tilting plateau
- Usable on two sides
- 4 wheels





### Wardrobe

- Dimensions 100 x 200 x 60 cm (W x H x D)
- Safe for valuables
- Lockable



# Configure your bed ... ... according to your wishes

## Step 1: Select the type of bed.

Scintilla 1000



X-scissor

Size: 90 x 200 cm

Height adjustment 38 – 80 cm

Scintilla 2000



Y-scissor

Size: 90 x 200 cm

Height adjustment 25 – 80 cm

## Step 2: Select the colour.

Piemont  
Mokka

Country  
Oak

Beech

## Step 3 + 4: Select the head- and footboard.

Scintilla 1000/2000/XL




Straight version

Scintilla 1000/2000



Arch version

## Step 5 + 6: Select the type of side rail.

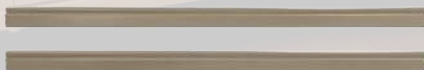

Full-length side rail






Separate side rail  $\frac{1}{2}$  +  $\frac{1}{2}$



Separate side rail  $\frac{1}{3}$  +  $\frac{1}{3}$



## Step 7: Select the type of remote control.

Comfort Care beds	Number of function buttons	Comfort setting	Anti-Trendelenburg-function	Trendelenburg-function
 Scintilla 1000	8	✓		
 Scintilla 2000	9	✓	✓	
 Scintilla 2000	10	✓	✓	✓

## Step 8: Select the type of wheels.



100 mm wheels



75 mm wheels



wheels (75 mm) + wheel covers  
(only available for the Scintilla 2000)

## Step 9: Select the brake system.

Standard: Axle braking system (2 by 2)

Optional: Central braking systems (4 wheels)

## Accessory



Bed side lamp



Automatically/manually adjustable trapeze bar



Comfort lying surface with panel system



Stand-up bracket



Underbed lighting

Further options:

- Battery
- Wall spacer
- Infusion stand

## Furniture

Various bedside tables



Various wardrobes





**Funke Medical GmbH**  
Ährenfeld 10  
46348 Raesfeld  
Tel.: +49-2865-266 93 50  
Fax: +49-2865-8133  
info@funke-medical.de  
www.funke-medical.de



Safety  
Quality  
Functionality



www.tuv.com  
ID: 0303003688